



MICROSCOPES



cs@irisanalytical.com



+1 706 760 4500



www.irisanalytical.com



Laboratory Metallurgical Microscope IMM-1000 Series is used to view samples that will fit on the microscope stage and it is the most common type of metallurgical microscope. This series comprises an infinite optical system and an illumination system.


Laboratory Metallurgical Microscope IMM-1100 Series


Laboratory Metallurgical Microscope IMM-1100 Series is a part of the range of laboratory microscope and stands out through its design which is ergonomic, vigorous and constant and stable. This laboratory microscope is used to view samples that will fit on the microscope stage. This range, with its large working distance with infinite plan achromatic: 4X/0.1/ ∞ -17.3 mm, 10X/0.25/ ∞ -10.0 mm, 40X/0.60/ ∞ /0-2.9 mm objective. Adjustable 6V/20W halogen lamp Koehler illumination and epi-illuminator ensures the optimum illumination of the materials to be tested. The illuminating system consists of a high-intensity light source, Abbe N.A. 1.25 condenser lenses with Iris Diaphragm. The binocular and trinocular viewing head with wild field WF10X18 mm eyepieces fitted with 48-75mm interpupillary distance as standard. Blue, Green, Yellow and Ground filters are often used to recover detail and polarizing filters, to produce surface glare and to recover grain boundary definition.



- Infinite plan achromatic objective, and Quadruple nosepiece
- Wild field S WF10X18 mm eyepieces
- Coaxial coarse & fine adjustment fine division of length 0.002 mm, coarse stroke 37.7 mm/rotation, fine stroke 0.2 mm/rotation, range of 20 mm
- Double layer mechanical stage of length 150×140 mm along with moving range 75×50mm
- Strain-free swing out condenser N.A.0.9/0.25
- 6V/20W Halogen lamp reflected and transmitted light source
- Different filters Blue, Green, Yellow and Ground filters
- Binocular 30° inclined and a 360° rotatable head with interpupillary distance 48 to 75 mm and trinocular 30° inclined and a 360° rotatable head with interpupillary distance 48 to 75 mm
- Small foot print and light weight
- Easy operating with ergonomic design
- Excellent image quality with Infinite optical system
- Optional accessories Specimen Presser For metallurgical specimen preparation; Polarizer and Analyzer; Photo and Video Adapter

Sample Image

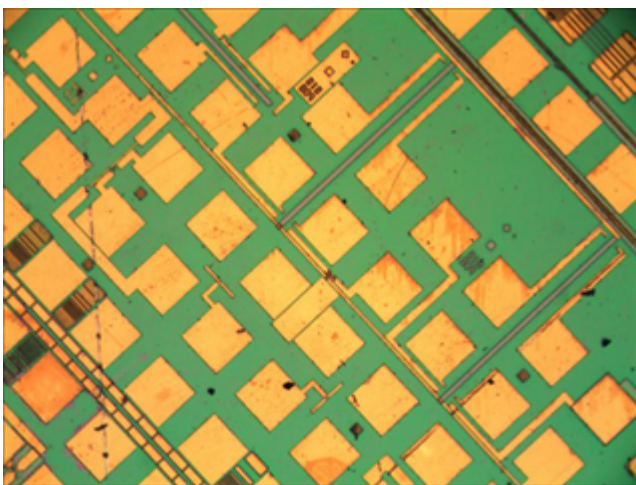

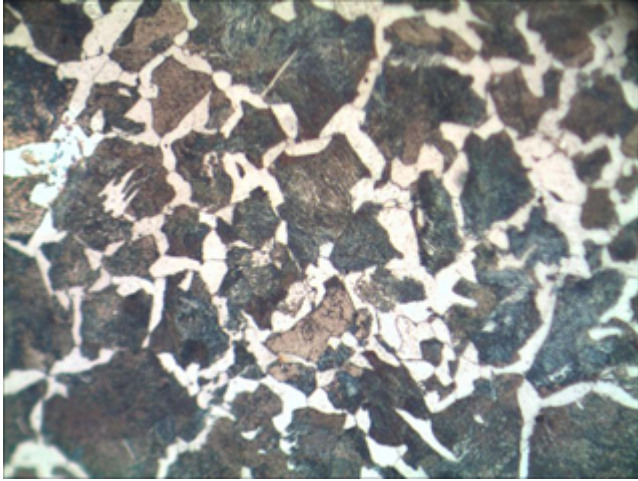

Laboratory Metallurgical Microscope IMM-1200 Series

Laboratory Metallurgical Microscope IMM-1200 Series is a part of the range of laboratory microscope and stands out through its design which is ergonomic, vigorous and constant and stable. This laboratory microscope is used to view samples that will fit on the microscope stage. This range, with its large working distance with infinite plan: 5X/0.10, 24-23mm, 10X/0.25-18.48mm, 20X/0.40 (S)-8.35mm, 50X/0.7-1.95mm objective. Adjustable 12V/50W reflected and transmitted halogen lamp to ensure the optimum illumination of the materials to be tested. The illuminating system consists of a high-intensity light source, Abbe N.A. 1.25 condenser lenses with Iris Diaphragm. The binocular viewing head with wild field WF10X22 mm, WF10X22 mm eyepieces fitted with 48-76 mm interpupillary distance as standard. Blue with optional Green, Yellow and Ground filters are often used to recover detail and polarizing filters, to produce surface glare and to recover grain boundary definition.



- Long Working Distance, Infinite Plan objective, and Backward Quintuple nosepiece
- Wild field Adjustable WF10X22mm, WF10X22mm eyepieces
- Coaxial Coarse and Fine Focusing, Fine Division 0.002mm, Fine Stroke 0.2mm per Rotation, Moving Range 25mm. Tightness adjustment, upper distance limiter
- Mechanical Stage: 200x150mm, Travel range: 77x52mm. Stage fast lowering unit:
- Abbe N.A. 1.25 with Iris Diaphragm
- 12V/50W Halogen lamp reflected and transmitted light source
- Different filters Blue with optional amber, green, grey filters
- Binocular 30° inclined and a 360° rotatable head with interpupillary distance 48 to 75 mm and trinocular 30° inclined and a 360° rotatable head with interpupillary distance 48 to 75 mm
- Designed with small foot print and light weight
- Easy operating with ergonomic design
- Excellent image quality with Infinite optical system
- Optional accessories C-mount Adapter, filters


Laboratory Metallurgical Microscope IMM-1300 Series


Laboratory Metallurgical Microscope IMM-1300 Series is a part of the range of laboratory microscope and stands out through its design which is ergonomic, vigorous and constant and stable. This laboratory microscope is used to view samples that will fit on the microscope stage. This range, with its large working distance with infinite plan achromatic: 5X/0.12/ ∞ /(BF&DF)-10 mm, 10X/0.25/ ∞ /(BF&DF)-10 mm objective. Adjustable 12V/50W reflected and transmitted halogen lamp to ensure the optimum illumination of the materials to be tested. The illuminating system consists of a high-intensity light source, Abbe N.A. 1.25 condenser lenses with Iris Diaphragm. The binocular viewing head with wild field WF10X22 mm, WF10X22 mm eyepieces fitted with 48-75mm interpupillary distance as standard. Green, Gray and ground glass filter are often used to recover detail and polarizing filters, to produce surface glare and to recover grain boundary definition.



- Long Working Distance, Infinite Plan objective
- Coaxial coarse & fine adjustment, fine division 2 μ m, moving range 30mm
- Offer stages that are designed to move small sample slides or entire microscopes including Mechanical Stage: to facilitates sample observations of length 200x150mm, Travel range: 77x52mm
- Wild field Adjustable WF10X22mm, WF10X22mm eyepieces
- 12V/50W Halogen lamp reflected and transmitted light source
- Different filters Green, Gray and ground glass filter
- Binocular 30° inclined and a 360° rotatable head with interpupillary distance 48 to 75 mm and trinocular 30° inclined and a 360° rotatable head with interpupillary distance 48 to 75 mm
- Small foot print and light weight
- Easy operating with ergonomic design
- Excellent image quality with Infinite optical system
- Optional accessories C-mount Adapter, filter

Sample Image

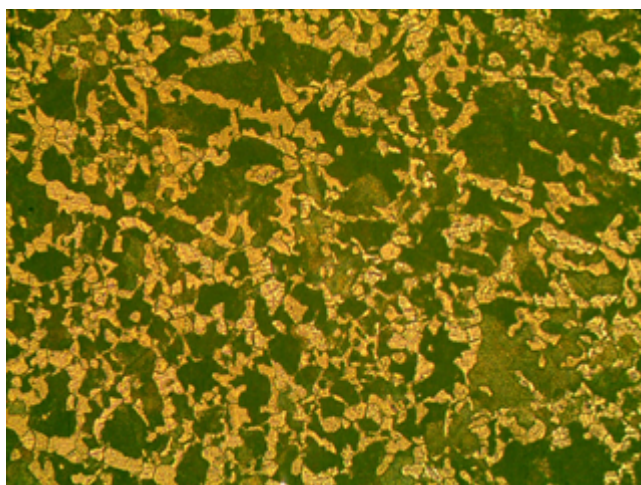

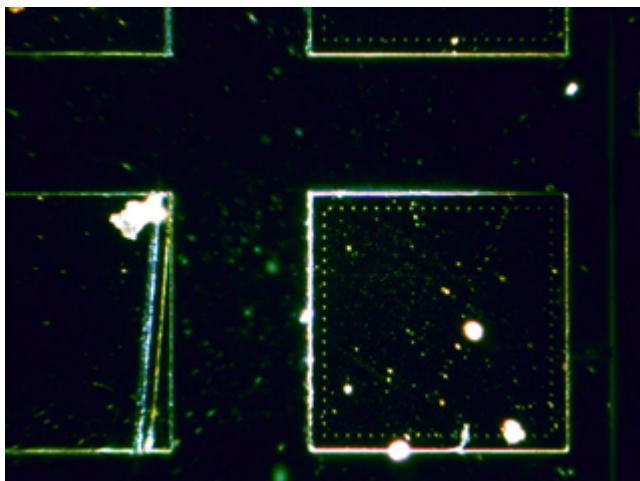

Accessories



