



MICROSCOPES



cs@irisanalytical.com



+1 706 760 4500



www.irisanalytical.com


Parallel Light Zoom Stereo Microscope IMS-3200 Series is the advanced Tilting type of microscope. The Tilting trinocular viewing head with point high field eyepieces and has excellent zoom objective to choose 0.8X-10X and Plan Apochromatic: 1X78 mm, with optional choice 2X20 mm objectives, and each eye tube has interpupillary adjustment. IMS-3200 Series is capable of providing strong images with a broad field of view, the multipurpose and perfect for basic studies of samples properties. Come with LED lights, allowing samples to be illuminated for a more comprehensive field of view. Suitable for Industrial, Medical and Scientific Research Areas, circuit board repair, circuit board inspection, surface mount technology work, electronics inspection, coin collecting, gemology and gemstone setting, engraving, and repair and inspection of small parts. Provide upright, un-reversed 3 D images. Used to examine the huge specimen under 0.8X-10X low zoom magnification objective and 1:12.5 zoom ratios.

- Plan Apochromatic: 1X78mm objective
- Wide Field Zoom Stereo Microscope
- Upright, un-reversed 3 D images
- High point wild field PL10x/23mm eyepieces
- Maximum magnification of up to 0.8X-10X
- Coarse and fine coaxial focus system, integrated body with focus holder, coarse range of length 50 mm, fine precision 0.002 mm
- 5W LED lamp illumination
- Binocular 45° inclined and a 360° rotatable viewing head with Interpupillary adjusting distance 55-75 mm and trinocular viewing head
- Small foot print and light weight
- Easy operating with ergonomic design
- Pillar stand
- Quality Image with infinite optical system
- Suitable for Industrial, Medical and Scientific Research Areas, circuit board repair, circuit board inspection, surface mount technology work, electronics inspection, coin collecting, gemology and gemstone setting, engraving, and repair and inspection of small parts
- LED ring light and camera adapter as an optional accessories



Accessories



Tilting viewing head for comfortable operation

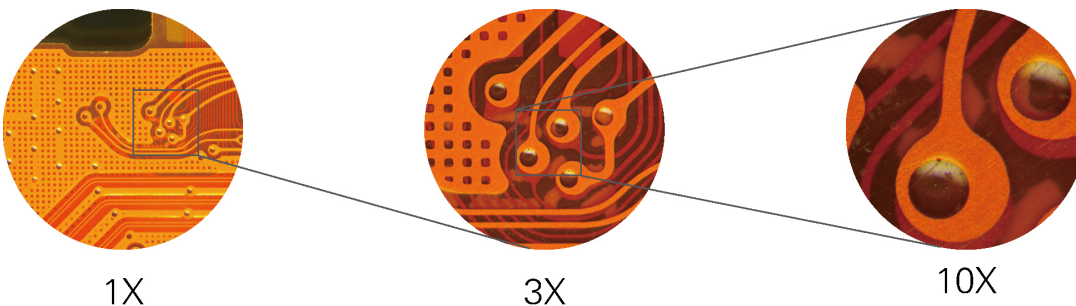

Large zoom ratio 12.5:1

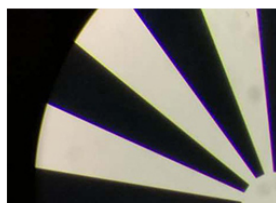
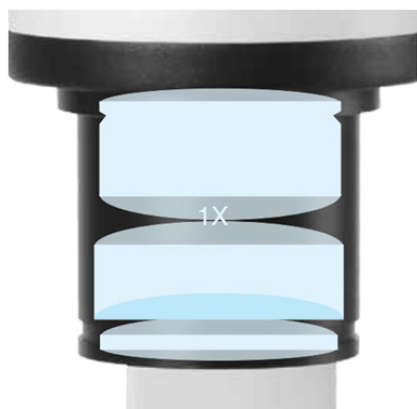


Other objective imaging



Apochromatic objective imaging

Apochromatic objective

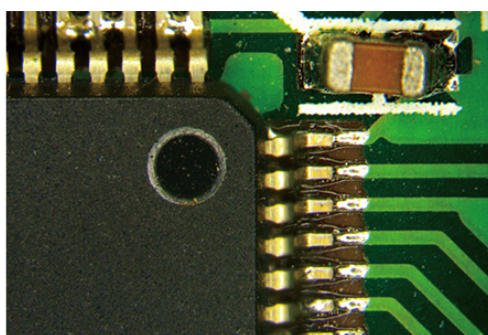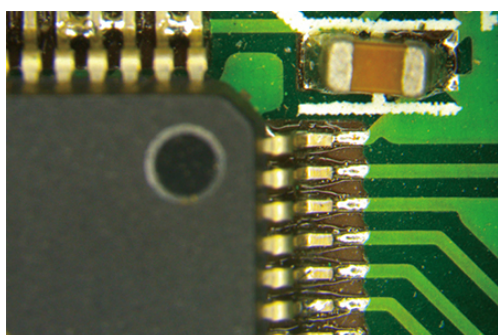

Aperture diaphragm

